

QuantideX® qPCR PML-RARA Kit*

Sensitive & Scalable Detection of Three Primary PML-RARA Fusion Transcripts

The QuantideX® qPCR PML-RARA Kit offers sensitive and scalable detection of three PML-RARA aberrant fusions in a single reaction, with results in less than 4 hours. Bcr1, bcr2, and bcr3 aberrant fusions of the PML-RARA gene are hallmarks of acute promyelocytic leukemia (APL), a distinct and aggressive subtype of acute myeloid leukemia (AML).

Reduced Complexity

- Ready-to-use design with all necessary RT-qPCR reagents, controls, and calibrators included
- Compatible with broadly installed laboratory equipment

Optimized Workflow

- Multiplexed design allows detection of three primary isoforms of PML-RARA (bcr1/bcr2/bcr3) and the reference gene in one reaction
- Automated data analysis and reporting software

Quality Results

- 4-point calibrator system spans assay dynamic range
- Highly sensitive limit of detection (~10⁻⁴)

FIGURE 1.

QuantideX qPCR PML-RARA Kit Workflow

	RT	qPCR	Data Analysis	Total
Hands-On Time	15 min	10 min	2 min	27 min
Instrument Time	1 hr 10 min	1 hr 30 min		<3 hrs



[Learn More](#)

Ordering Information

Part #
A00598

Product
QuantideX qPCR PML-RARA Kit


Reactions
60

*For Research Use Only. Not for use in diagnostic procedures.
Trademarks and registered trademarks are the property of their respective owners.

3000-371 V3_0925

biotechne® / Global Developer, Manufacturer, and Supplier of High-Quality Reagents, Analytical Instruments, and Precision Diagnostics.

INCLUDES R&D Systems™ Novus Biologicals™ Tocris Bioscience™ ProteinSimple™ Asuragen™ ACD™ Lunaphore™

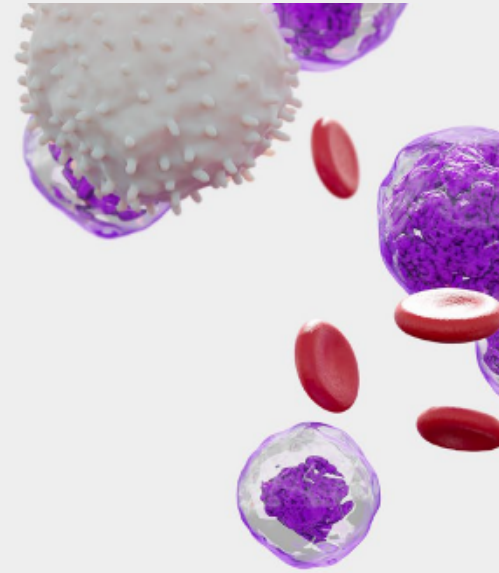

Distributed by Abacus dx

1800 ABACUS (AUS) 0800 222 170 (NZ) | info@abacusdx.com | www.abacusdx.com

abacus dx