


Lab & Production Materials

Millipore Sigma


Cell Marque™ Tissue Diagnostics

ATRX Rabbit Polyclonal Antibody For use in neuropathology applications


Diffuse gliomas are classified based on histological and molecular features to achieve an integrated diagnosis. Molecular diagnostic markers include IDH mutation, 1p/19g co-deletion, and TP53 mutation. ATRX is a chromatin remodeling protein, and its mutation status may be used as a molecular diagnostic marker within the diffuse glioma classification algorithm. Anti-ATRX is used to identify mutant ATRX by a loss of ATRX expression in neoplastic cells when compared with internal positive controls (endothelial cells, glia, and neurons). Grade II/III astrocytoma classification includes IDH mutant, ATRX mutant, and 1p/19g retention, while grade II/III oligodendroglioma includes IDH mutant, ATRX wildtype, and 1p/19q co-deletion; p53 expression may also serve as an aid in diagnosis. ATRX mutation is frequently, but not always, mutually exclusive with 1p/19q co-deletion.¹⁻⁶

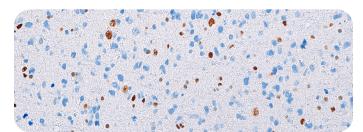

Astrocytoma

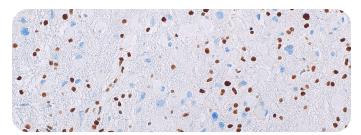

Astrocytoma

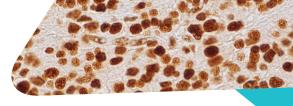

Oligodendroglioma

Ordering Information:

Description	Cat No.
0.1 mL concentrate	485A-14
0.5 mL concentrate	485A-15
1.0 mL concentrate	485A-16
1.0 mL predilute ready-to-use	485A-17
7.0 mL predilute ready-to-use	485A-18



Intended Use:

This antibody is intended for *in vitro* diagnostic (IVD) use. ATRX Rabbit Polyclonal Antibody is intended for laboratory use in the detection of the ATRX protein in formalin-fixed, paraffin-embedded human tissue stained in qualitative immunohistochemistry (IHC) testing. The results using this product should be interpreted by a qualified pathologist in conjunction with the patient's relevant clinical history, other diagnostic tests and proper controls.

Product Information:

Visualization: Nuclear

Controls: Astrocytoma, Glioblastoma, Brain, Oligodendroglioma Dilution: 1:50

Associated Specialty: Neuropathology

References:

- 1. Cancer Genome Atlas Research Network. Comprehensive, Integrative Genomic Analysis of Diffuse Lower-Grade Gliomas. N Engl J Med. 2015 Jun 25;372(26):2481-98.
- Reuss DE, et al. ATRX and IDH1-R132H immunohistochemistry with subsequent copy number analysis and IDH sequencing as a basis for an "integrated" diagnostic approach for adult astrocytoma, oligodendroglioma and glioblastoma. Acta Neuropathol. 2015 Jan;129(1):133-46.
- Louis DN, et al. The 2016 World Health Organization Classifi cation of Tumors of the Central Nervous System: a summary. Acta Neuropathol. 2016 Jun;131(6):803-20.
- Ikemura M, et al. Utility of ATRX immunohistochemistry in diagnosis of adult diff use gliomas. Histopathology. 2016 Aug;69(2):260-7.
- Yamamichi A, et al. Immunohistochemical ATRX expression is not a surrogate for 1p19q codeletion. Brain Tumor Pathol. 2018 Apr;35(2):106-113.

MilliporeSigma 400 Summit Drive

Burlington, MA 01803

www.SigmaAldrich.com

 Wood MD, et al. Applications of molecular neuro-oncology - a review of diff use glioma integrated diagnosis and emerging molecular entities. Diagn Pathol. 2019 Apr 9;14(1):29.

The product is FDA Class 1 in the US. The product featured is not available in all countries.

Contact your local sales representative or distributor for details.

Toll Free: 800.665.7284 Phone: 916.746.8900 Fax: 916.746.8900 Email: service@ccellmarque.com www.cellmarque.com

Copyright © 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. MilliporeSigma, the vibrant M, Sigma-Aldrich and Cell Marque are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

Lit. No. MS FL8635EN Ver. 1.1 38559 05/2022

Distributed by Abacus dx

1800 ABACUS (AUS) 0800 222 170 (NZ) | info@abacusdx.com | www.abacusdx.com

