
# PermaLife<sup>™</sup> Cell Culture Bags

Cell Culture and Cell Expansion

Cryopreservation

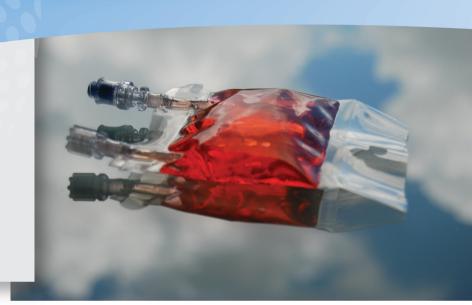

#### Advantages of the OriGen PermaLife™ Cell Culture Bags

The OriGen PermaLife<sup>™</sup> Cell Culture bags represent an optimal way to protect, store, and freeze your critical tissue and cell cultures. Made of inert fluoro-ethyl polymer (FEP), the PermaLife<sup>™</sup> bag contains no plasticizers or chemical residuals which could leach into and alter delicate cell cultures or medium.

PermaLife<sup>™</sup> FEP bags are a closed system, so the possibility of contamination is greatly reduced. Swabbable, needle-free valves are provided on every bag to preserve sterile access and reduce the chance of needle sticks. Most PermaLife<sup>™</sup> bags are sold with one female luer and one needle-free valve. The PL07 only has one port with one needle-free valve.

PermaLife<sup>™</sup> bags are gas permeable and provide a high rate of oxygen, carbon dioxide, and nitrogen transfer to cells under culture. Additionally, water permeability is very low; water loss is negligible and humidified incubators – potential sources of mold and contamination - are not required.

The PermaLife<sup>™</sup> bag has an optical transmission of at least 96%. Observation of the cultures under a microscope is possible without risking contamination of the culture.


The PermaLife bag can be utilized for cryogenic preservation.




## **PermaLife**

### **Cell Culture Bags**

#### **Technical Information**

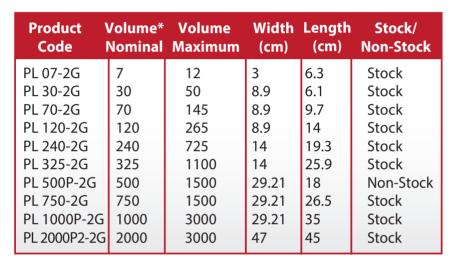

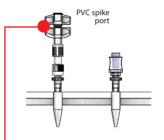

#### **Table Legend:**

"P" in the catalog number denotes it includes a PVC spike port.

Length & Width of active culture compartment only.

\* Nominal Volume at 1cm fill thickness, volume in ml



#### **Special Orders Available**

Product configurations which include a modification of labels, different connectors and cap materials, multi-port connections, SCD tube sets and more are available. Let OriGen know your needs and how we can help improve productivity. Minimum quantities apply.

#### **Ordering Information:**

Sterile, for single use only. Packed and sold individually. To order, contact the office near you. For availability in your country, contact your OriGen representative.



